

set[®] feeding





set[®] feeding

Enteral nutrition requires equipment that is dependable in every routine. Our feeding portfolio supports safe and consistent enteral delivery across clinical settings.

With products designed for routine use in hospitals, pharmacies, and home-care environments, this category provides a clear overview of our feeding syringes and related accessories.



We care.

Feeding Instruction


What Is a Feeding Syringe?

ENFit® Tip

The ENFit® connector is a standardized connection designed for enteral feeding systems that improves patient safety by preventing misconnections.




The purple plunger that enables differentiation.

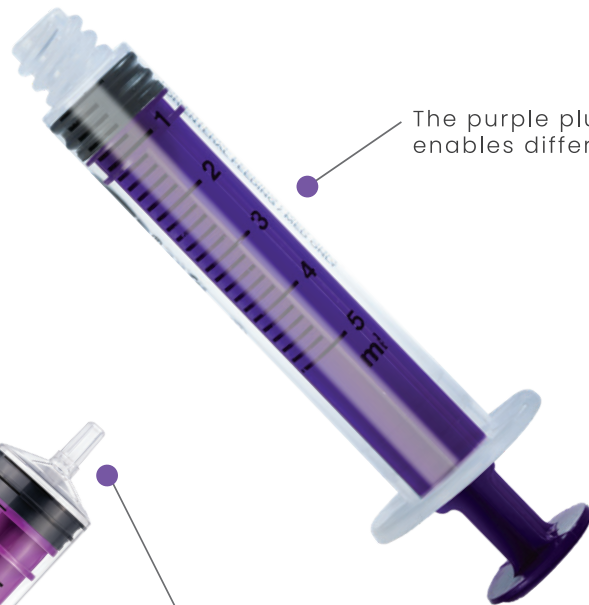

Oral Tip

Oral tip syringe is designed for oral feeding use, helping prevent misconnections and ensuring safe, accurate administration.

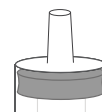





Set® Oral

Features

- Designed for safe dispensing and oral administration of liquids, featuring a self-righting cap design
- Smooth plunger movement with silicone lubrication
- Latex-free gasket ensuring reliable sealing
- Clear and accurate graduation for precise dosing
- Available in 1 ml, 3 ml, and 5 ml sizes ensuring enhanced micro-dosing accuracy
- Transparent barrel offering easy visibility of contents
- Single-use device ensuring patient and user safety

Packaging

- Sterile single blister pack
- Non-sterile options
- Without or with cap.


Customization

- Plunger color options 
- Clear or Amber barrel options 

Conformance & Compliance

- ISO 7886-1 – ISO 80369-1
- EU 2017/745 (MDR)
- EO sterilized
- Shelf life: 5 years

Available Sizes & Types

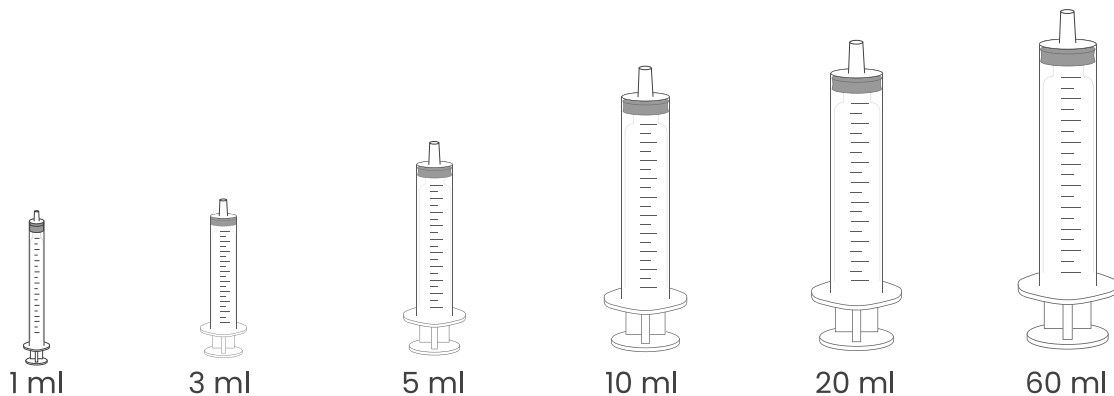

Product Table

Size	Graduated Capacity	Inner Box	Carton Qty
1 ml	1 ml	100/500 pcs	3200/3000 pcs
3 ml	3 ml	100/350 pcs	2400/2800 pcs
5 ml	5 ml	100/250 pcs	1600/2000 pcs
10 ml	10 ml	100/150 pcs	1200/1200 pcs
20 ml	20 ml	100 pcs	600 pcs
60 ml	60 ml	60 pcs	240 pcs

Set® ENFit®

Features

- Designed for precise enteral medication delivery and fluid measurement
- Smooth and controlled plunger movement with silicone lubrication
- Latex-free gasket ensuring excellent sealing
- Clear, easy-to-read graduation for accurate dosing
- Available in 1 ml, 3 ml, and 5 ml sizes ensuring enhanced micro-dosing accuracy
- Fully compatible with ENFit® Standard and ENFit® Low Dose connectors
- Transparent barrel for clear visibility during feeding and dosing
- Single-use device ensuring patient and caregiver safety

Packaging

- Sterile single blister pack
- Non-sterile options
- Without or with cap.


Customization

- Plunger color options 
- Clear or Amber barrel options  

Conformance & Compliance

- ISO 7886-1, ISO 7886-2 - ISO 20695 - ISO 80369-1 - ISO 80369-3
- EU 2017/745 (MDR)
- EO sterilized
- Shelf life: 5 years

Available Sizes & Types

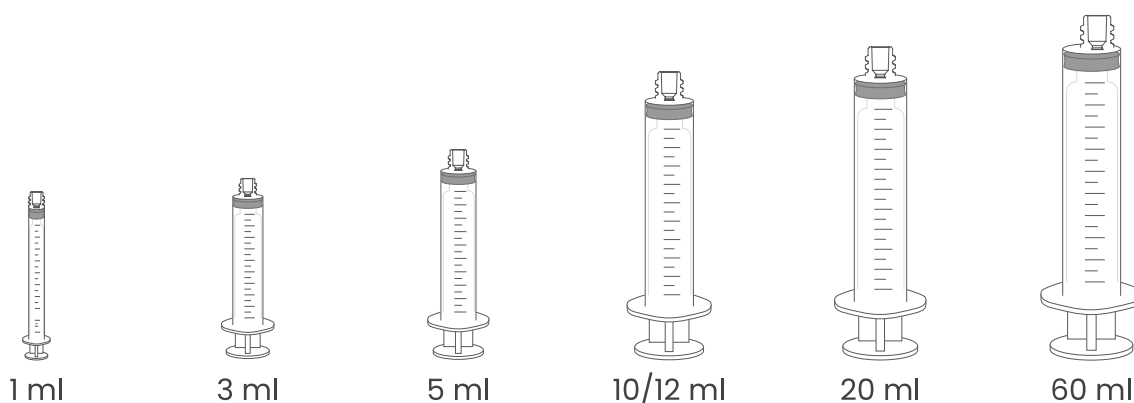

Product Table

Size	Concentric	Eccentric	Low Dose Tip	Standard Tip	Graduated Capacity	Inner Box	Carton Qty
1 ml	•		•		1ml	100/500 pcs	3200/3000 pcs
3 ml	•		•		3ml	100/350 pcs	2400/2800 pcs
5 ml	•		•	•	5ml	100/250 pcs	1600/2000 pcs
10/12 ml	•			•	10ml	100/150 pcs	1200/1200 pcs
20 ml	•			•	20ml	100 pcs	600 pcs
60 ml		•		•	60ml	60 pcs	240 pcs